



ampacdental

DENTAL EQUIPMENT SPECIALISTS

2022 Catalogue



swident



NSK
Create it





LED light

Assistant table

- Luzzani or DCI triplex syringe
- HVE and LVE hoses
- Touch panel for the operation of cup filler, bowl flush, treatment positions and light activation.

Ceramic spittoon

removable and autoclavable spittoon for easy cleaning. Rotates 180° for convenient patient access.

Patient chair

- Multi-functional foot pedal for easy control of treatment positions
- Emergency cut off switch
- Easy accessible suction filter
- Trendelenburg position
- Chrome or white armrests
- Large selection of vibrant colours
- Optional clean water system



"The Swident chairs have been absolutely incredible"

MADE IN ITALY


- Dr Kanowitz
Principal Dental Surgeon, Hornsby

Doctor table

- 2 fibre optic high speed lines
- 1 electric micro motor
- optional scaler
- choose from the following:



traditional or continental

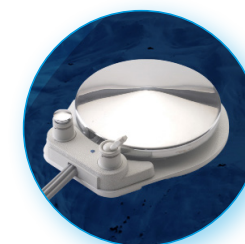

pneumatic



electric



DCI



Luzzani

Friendly to use - Easy on wallet

- Built with the budget conscious professional in mind, the Swident **Friend Easy** boasts a simple design in a robust, yet compact frame that fits into operatories with limited space. Includes Faro LED Maya dental light.
- Pneumatic foot control allows progressive regulation of the micro motor and Spray and Chip blower function.
- **Traditional or Continental delivery system.** Equipped with: 1 high speed fibre optic line, 1 low speed line and space for scaler. Select from either DCI or Luzzani triplex syringe.
- Available as an orthodontic unit without the spittoon and doctor table.



"The Friend Easy fits perfectly in our modest surgery room that measures only 3 x 3 metres, still affording extra room for the operator and assistant"

- Dr Praveen Ayethepalle, Perth Dentist

Now available in Rear Suction

- Enjoy all the benefits of the Swident Partner chair, now with an optional rear suction functionality.
- The pneumatic dental chair includes Faro Maya LED dental light, pneumatic foot pedal, 2 high speed lines, 1 low speed line and optional scaler.
- Rear suction arm includes triplex syringe, HVE, LVE.
- Soft memory foam upholstery in various vibrant colours.



MADE IN ITALY



New

Swident cart



Cart inclusions

- Luzzani or DCI triplex syringe
- 1 highspeed FO, 1 lowspeed line
- 2 additional spaces
- Pneumatic foot control
- Optional clean water system

Doctor and assistant stools



DCI Delivery Systems

Suitable for orthodontic units, left and right hand use or surgeries with limited space

DCI delivery systems include:

- DCI Precision Comfort triplex syringe
- 3 handpiece lines (eg. highspeed, lowspeed line, space for scaler)
- Individual drive air and water coolant controls
- Wet/Dry foot control



Mobile cart system #4511

- Sturdy all metal control head with smooth lines for easy cleaning
- Stainless steel tray and holder
- Deluxe clean water system
- 1/4" water QD w/flow control
- 7ft of 1/4" umbilical
- Optional fibre optic lines



Unit or cabinet mount #4478

- Available in left or right hand configuration
- Flex arm system mounts on 2" post
- Deluxe clean water system
- 7ft of 5/8" umbilical
- Fitted to A-dec® and other equipment brands
- Optional fibre optic lines
- Optional storage area for high tech accessories



DCI Delivery systems



DCI III Manual Unit
#4502




Panel mount automatic control for
3 handpieces #4416



Cabinet mount #4370
Unit mount #4376


DCI Repair and Service Kits

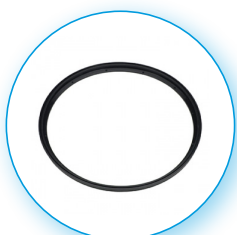

Gasket, to fit A-dec® W&H
Lisa MB17 sterilizer #2186



1/4" poly Q.D panel mount
with shut-off, female #0011

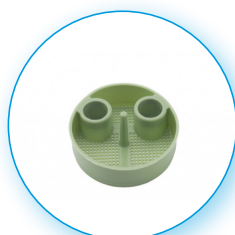

Solids collector disposable
screen 1-7/8" diam. #5818



Service kit, to fit A-dec® Century
Plus® control block #9154



Quick clean syringe buttons
and repair kit #3098



O-ring Pk12 #2203
corresponds with
A-dec® 030-006-00


DCI Dental Chair components



A-dec® valve #5062



Water bottle 2l #8164
1l #8669 750ml #8128

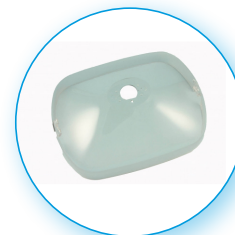

Light shield for A-dec®
halogen lights #9390



Belmont® Clesta light shield
#8959




Foot control wet/dry
& chip blower #6315



Kavo® LED globe #9372



Belmont Clesta® globe #8890

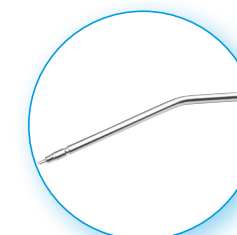

Syringe tip #3060 autoclavable
Syringe tip A-dec® #3055



Syringe, one button #3640
no tubing or kit, non autoclavable



Syringe, Precision Comfort #3604
non autoclavable



Syringe, autoclavable #3370
valve core, head only

These parts are manufactured by DCI to fit A-dec® equipment. DCI has no affiliation with A-dec Inc.
DCI replacement parts are not manufactured under license or supervision of A-dec®

Nardi compressors are manufactured to the highest standard to deliver quality, OIL FREE, DRY, CLEAN AIR



25L 73 dB



25L hood 67 dB



30L 64 dB; 50L 62 dB



- Oil Free, Dry, Pure Air
- Maintenance Free
- 3 year parts warranty
- Made in Italy
- Certified by Workcover NSW
- Suitable for HOT climates



50L hood 67dB



50L 73 dB

Oil residue can damage handpieces and fine pneumatic controls in the dental units. An oil free dental compressor is vital in protecting the health of your staff and patients. It will ultimately increase the longevity of your equipment.

Tanks are furnace painted on the inside and outside using epoxy electrostatic dust to prevent corrosion. All the compressors are certified by Workcover NSW.



50L 73 dB



90L 73 dB



90L CAD/CAM 73 dB

| tank size | dryer | cylinders | air intake l/min | output l/min @5bar | pressure Bar | power HP | decibels dB | weight Kg | dimension LxWxH cm |
|--------------|-------|-----------|---------------------|--------------------------|-----------------|-------------|----------------|--------------|-----------------------|
| 25 l | yes | 2 | 350 | 156 | 8 | 2.00 | 73 | 44 | 59x43x66 |
| 25 l hood | yes | 2 | 350 | 156 | 8 | 2.00 | 67 | 48 | 70x43x78 |
| 30 l cabinet | yes | 2 | 350 | 156 | 8 | 2.00 | 64 | 79 | 71x44x76 |
| 50 l W | yes | 2 | 350 | 156 | 8 | 2.00 | 73 | 51 | 84x44x68 |
| 50 l T | yes | 4 | 480 | 240 | 8 | 2.50 | 73 | 60 | 64x45x89 |
| 50 l hood | yes | 2 | 350 | 156 | 8 | 2.00 | 67 | 55 | 70x43x98 |
| 50 l cabinet | yes | 2 | 350 | 156 | 8 | 2.00 | 62 | 96 | 85x50x82 |
| 90 l | yes | 4 | 480 | 240 | 8 | 2.50 | 73 | 86 | 110x59x81 |
| 90 l cadcam | 2 | 4 | 480 | 240* | 8 | 2.50 | 73 | 86 | 110x59x81 |

*output @ 8bar=190 l/min

MADE IN ITALY



Intra-oral X RAY



State-of-the-art technology, superior performance and reliable service makes the Trident **RiX-70 DC** the clear choice.

Leading edge .04mm focal spot technology delivers outstanding image quality. Whether your using film or digital applications, your images are sharper with incredible detail and superior contrast. The result is improved diagnostics and enhanced patient care.

The main control panel has easy to read icons offering numerous preset exposure settings for patient size (child, small adult, large adult), tooth structures, film type (plates, sensor or film). An additional keypad is installed near the tube head to manually adjust the exposure time and tooth type.

Adjustable tube voltage of 60-70 kV combined with steady current of 7mA allows the operator to select various settings to produce detailed images for improved sharpness and contrast.

With a strong, robust design, the wall mounted xray arm is available in lengths of 100cm, 80cm and 60cm to suit even the most challenging of spaces.

Backed up by 2 year warranty on parts and 1 year on labour, the Rix-70 DC xray represents excellent value for money. All products are engineered and manufactured by Italian company Trident, who specialise in radiological products.



RiX-70 DC

OPG Panoramic 2D



Provide superior treatment with the **Panoramico 2D** OPG. An evolutionary machine with a compact streamlined design that fits easily into your practice.

Achieve outstanding image quality when producing any of the following programs:

- Adult and child panoramic
- Bitewing
- Sectoral panoramic and reduced radiation panoramic
- Orthopantomogram
- TMJ
- Paranasal sinuses



OPG Panoramico 2D

Phosphor Plate system



Capture, scan and share clear digital images with the **Reader** Digital Phosphor Plate system.

Fast readout time 10 seconds. Plates are cleaned and ready for immediate use.

Support up to 10 workstations. TWAIN driver available. Integrates with most popular practice management software.

Starter pack includes phosphor plate sizes 0,1,2,3, covers and sleeves.



Reader

Ultrasonic Cleaners

Warranty 2½ years circuit board, 5 years transducers, 10 years tank



L&R SweepZone AgION

- **Antimicrobial-coated housing** provides a continuous shield from bacteria, mould and fungi to keep your equipment clean and safe.
- **SweepZone technology** sweeps through the entire tank resulting in a total clean of all instruments irrespective of their position in the tank. Eliminates hot spots and weak areas!
- **Digital touch pad** allows easy control of the functions.
- **Stainless steel housing and lid** means it is corrosion, rust and fingerprint resistant.
- **4L** LxWxH 32.2 x 17.5 x 22.2; 67dB
12.3L 41.9 x 25.4 x 30.5; 67dB



4L SweepZone



Recessed model

Benchtop model

- Flush lid design provides a clean streamline benchtop area.
- Recessed cover collects condensation.
- Recessed model sizes: 12.3L, 13.6L, 24.6L

New

Ultradose Wavecheck ultrasonic test monitor



Strips measure these variables:

- Cavitation
- Time and Temperature
- Solution ratio

Wavecheck® is a ready to use indicator, specially designed to test the efficacy of your ultrasonic cleaner.

Steam Steriliser

Warranty 2 years or 2500 cycles



The Melag brand is synonymous with quality, economy, speed and above all unparalleled sterilisation treatment for your dental instruments. Proudly made in Germany, Melag steam sterilisers are manufactured and tested under the strictest conditions to ensure complete customer satisfaction.

Evolution B Class autoclave incorporates revolutionary patented technologies to deliver faster, secure and more efficient processing of your instruments.

Double Jacket technology ensures air is rapidly removed from all load types (solid, hollow, wrapped and unwrapped) in the sterilisation chamber. The result is fast and efficient processing times.

Innovative and patented DRYtelligence reduces drying times by 80%. A special sensor automatically adapts to the drying load. A smaller load requires less drying time therefore limiting the exposure of the instruments to unnecessary heat and damage.

World's largest colour touch display is now even bigger. The intuitive menu structure enables program and cycle selection. Sterilisation run, log and other information is displayed in a clear, easy to read fashion.

Plumbed steriliser 22L B Class

- **20 minute Quick program B cycle** for a maximum of 1.5kg wrapped or 7kg unwrapped instruments.
- **30 minute Universal program** up to 7kg including drying.
- **Start time preselection** allows you to preselect and run a program at a later determined time.
- **Water treatment Meladem 40** ion exchanger unit automatically fills the steriliser with demineralised water and disposes of unused water.



Melag 44B+



Thanks to the increased chamber size you can select up to 8 trays and various configurations to suit your individual requirements.



If you want to label and trace the instruments following the successful sterilisation process, use **MELAp rint 60**.

Robust and Reliable. Patented technology to Guarantee 100% Air and Water separation



Teravac 6000
4 Dental Units

DMEGA suction motors use a Semi-Wet System. Both suction and separation are performed within the suction motor. The result is a significant reduction in costs as there is no need to install a water separator in the dental unit.

DMEGA Teravac 6000 is the most powerful and responsive suction motor with variable speed, high airflow and amalgam retention unit that enables up to 4 operators to work simultaneously. The intelligent Mitsubishi inverter monitors the vacuum levels and adapts accordingly to provide maximum suction whether you are running 1 or 4 chairs.

The patented innovative design uses 2 cyclone impellers to guarantee 100% air and liquid separation. Furthermore, the implementation of the patented Triple Wing System and auto drain valve, releases residue from the separator thereby reducing blockages, maintenance and ultimately costly repairs.

An optional amalgam retention device, available for all motors.



S101
1 Dental Unit

| | S101 | Teravac 3000 | Teravac 5000 | Teravac 6000 |
|--------------------|---------------|--------------|--------------|--------------|
| Treatment units | 1 unit | 2 units | 3 units | 4 units |
| Frequency | 50-60 Hz | 0-70 Hz | 0-70 Hz | 0-65 Hz |
| Main voltage | 230 VAC | 230 VAC | 230 VAC | 230 VAC |
| Max Airflow | 280 l/m | 600 l/m | 900 l/m | 1200 l/m |
| Max Vacuum | 150-180 m bar | 370 m bar | 370 m bar | 370 m bar |
| Decibels | 55 dB | 59 dB | 61 dB | 61 dB |
| Weight | 20 kg | 32 kg | 33 kg | 37 kg |
| Dimension cm WxLxH | 33x54x62 | 41x51x61 | 41x51x61 | 41x51x63 |

Knee break dental chair - Linda Next

Linda Next knee break chair provides ease of use and outstanding comfort for both the patient and doctor.

- Extremely easy to use thanks to the joystick foot control and built in programs for patient positioning, rinse position, last position, Trendelenberg and zero positions
- 3 operator profiles with 3 programmable positions
- Minimum height 500mm and maximum height of 875mm allows for standing procedures
- Synchronised chair movements allow smooth operation of the backrest and footrest. Slow mode feature option allows for precise patient positioning
- Articulated headrest accommodates different procedures
- Emergency safety switch and patient lifting capacity of 300Kg



Options:

- Rotating base; wheels
- Swiveling left and right armrests
- Soft upholstery, child booster seat
- Different headrests
- Backrest motorised height adjustment



Optional child booster seat

MADE IN ITALY



Small equipment



Alginate mixing machine



Work with perfectly mixed, bubble-free alginate every time. Brilliant durable little benchtop machine. Designed for use with any alginate or stone, this robust simple device mixes three scoops in just 11 seconds with no dust or waste. Never hand mix again.

"I bought 2 of these machines for my practice.. I was so impressed I now own 14!"

- Dr Joseph Geenty
Wollongong Orthodontist

Cordless LED Curing Lights

BlueLEX LD-M4

- Full power mode reaches 1000 mW/cm².
- Charge with micro usb, Wavelength 420-490nm
- Long battery life



BlueLEX GT-1200

- 3 light output modes: full power, ramp up, pulse
- Powerful output of over 1200 mW/cm²
- Built in radiometer



Capsule Mixer



a.max AM-1

- Stable, quiet and vibration free
- Powerful speed of 4800 opm ensuring consistent and homogenous mixing results
- Compatible with all capsules on the market including amalgam, glass ionomer, cements
- 2 memory settings, adjustable timer 1-30 sec.

Hoses



EMS® EMS® no pain
EMS® for Anthos®



Satelec® scaler hose



Fibre optic, non fibre optic hose



Kavo® fibre optic hose




Sirona® high speed hose




Sirona® micromotor hose

LVE

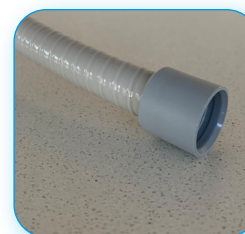

Cattani® 11mm
OEM 8mm, 10mm

HVE

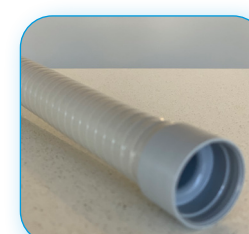

Cattani® 17mm
OEM 15mm, 16mm, 20mm

Connectors



Cattani® connectors and
terminals for HVE and LVE



Mini air micro motors
S-Max M205



Mini electric micro motors
NLX PLUS/MULTIPAD



Electric micromotor system
NLZ E



Couplings
PTL-CL-LED



Contra angle handpieces
Ti-Max Z95L SWITCH



High speed handpieces
S-Max M900L



Oral hygiene
Varios 170



Maintenance
iCare C2



Ask us about
Handpiece Repairs

Showrooms

Sydney, Perth, Newcastle, Townsville



ampac dental
DENTAL EQUIPMENT SPECIALISTS

30 YEARS



Free Call 1800 061 343 • T 02 9567 3555 • F 02 9599 2153
Suites 5&6, 39-45 George Street, PO Box 554, Rockdale NSW 2216
www.ampacdental.com.au • info@ampac-dental.com.au