

NARDI **COMPRESSORI** **MEDICAL**


Quality system certified
to ISO 9001:2015



NEW COLOR
RAL 7016
COMING SOON

MEDICAL COMPRESSOR
OIL-LESS

Catalogue

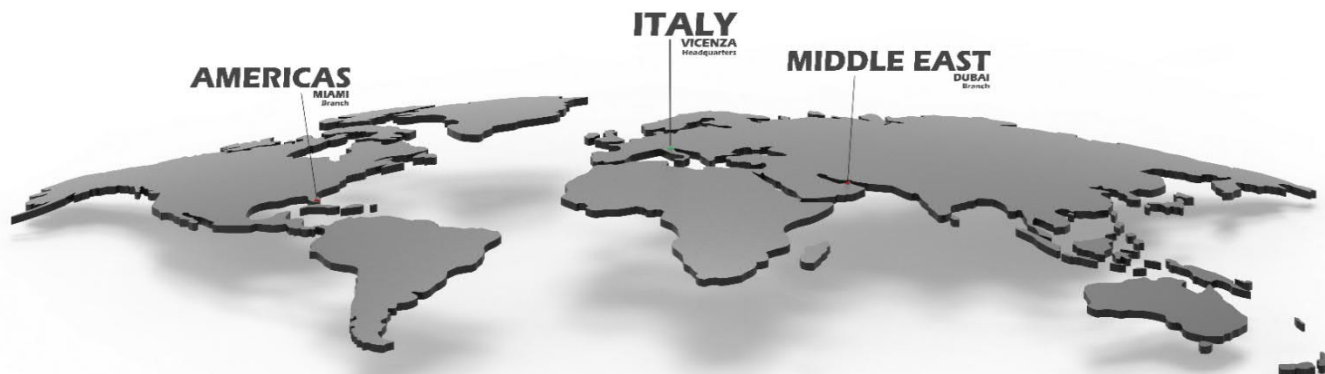

(ENG) NARDI COMPRESSORI is an Italian company that has focused on the design and manufacture of piston compressors and boosters since its very beginning.  The logo shows the two surname initials of the company founders: Nardi and Ghiotto. Over the years, the company has experienced exponential growth in the compressed air sector, which has led to the opening of two branches, located in Dubai (UAE) and Miami (USA). Today NARDI counts over 50 employees and collaborators and is appreciated in more than a hundred countries around the world thanks to its extensive distribution network. The flagship of this manufacturing company is the ability to study, design and project their products in-house, thanks to the technical expertise of CEO, Mr. Marco Ghiotto, and his technical staff. Production is entirely carried out in Italy, at the production facility in Montecchio Maggiore, and each unit is fully tested before being placed on the market, making the NARDI brand synonymous with quality. This commitment is proven by being UNI EN ISO 9001:2015 certified by TÜV. NARDI offers a wide and heterogeneous range of compressors, with models covering working pressures from 8 up to 500 bar, boasting a production of over 60,000 units per year.

(ITA) NARDI COMPRESSORI è una realtà italiana focalizzata sin dalla sua nascita, sulla progettazione e produzione di compressori e booster a piston.  Il logo rimanda alle iniziali dei cognomi dei fondatori dell'azienda: Nardi e Ghiotto. Negli anni l'azienda ha vissuto una crescita esponenziale nel settore dell'aria compressa, che ha permesso l'apertura di due filiali, site a Dubai e Miami. Ad oggi conta oltre 50 dipendenti e collaboratori ed è apprezzata in più di cento Paesi in tutto il mondo grazie alla sua vasta rete di distribuzione. Il fiore all'occhiello di questa azienda manifatturiera è la capacità di riuscire a studiare, disegnare e progettare i propri prodotti internamente, grazie alla preparazione tecnica del CEO, il Sig. Marco Ghiotto, e del suo staff tecnico. La produzione è totalmente realizzata in Italia, presso lo stabilimento di Montecchio maggiore, ed ogni unità viene completamente testata prima di essere immessa nel mercato, rendendo il marchio NARDI sinonimo di qualità. Tale impegno è comprovato dal fatto di essere certificata UNI EN ISO 9001:2015 dall'ente TÜV. NARDI offre una gamma vasta ed eterogenea di prodotti, con modelli che coprono pressioni di lavoro dagli 8 fino ai 500 bar, vantando una produzione di oltre 60.000 unità l'anno.



(ENG) The NARDI COMPRESSORI MEDICAL full range is designed and developed by expert technicians specifically for the dental and medical industry. The air tanks of this whole compressor line are painted internally with corrosion inhibiting epoxy electrostatic powder and the compressors are equipped with dryers that absorb humidity and produce a moisture-free and clean air quality. NARDI also avails of antibacterial filters to meet international standards required by the niche markets we cater for. The materials and components used to manufacture these modern, strong and silent compressors guarantee great performance and reliability and they are all certified as CLASS IIA MEDICAL devices and recognised by the Italian Ministry of Health.



MADE IN VENETO



EXCELLENCE & HIGH TECHNOLOGY

(ITA) La linea MEDICAL di NARDI COMPRESSORI è una gamma progettata e studiata appositamente per soddisfare gli standard richiesti dal settore dentale e medicale, offrendo compressori moderni e resistenti adatti ad essere utilizzati all'interno di laboratori odontotecnici, cliniche e studi dentistici. La particolarità di questi modelli, disponibili con diversi motori e telai, sta nel serbatoio, che viene sottoposto ad una verniciatura speciale che ne evita la corrosione interna nel tempo, oltre alla possibilità di installare accessori come essiccatori ad assorbimento o filtri antibatterici. Tutti i modelli della linea MEDICAL sono certificati CLASSE IIA e registrati come dispositivi medici presso il Ministero della Salute.



Packing Imballo

Weight with box (ESPRIT 60) : 26 Kg - 57,3 Lbs
Box dimensions cm : 49 x 49 x 72 (0,017m³)
Quantity per PALLET (110x110) : 8 pc.



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
ESPRIT 60/4	1	X	X
ESPRIT 60/2	1	X	X

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Revolution	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Giri	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	Rpm	dB	Kg - Lbs	cm (LxWxH)
ESPRIT 2V 60/4 25L	25	X	2	105 - 3,7	40 - 1,4 (5 BAR) 20 - 0,7 (8 BAR)	8 - 114	0,75 - 0,5	Single	1400	52	23 - 51	41 x 45 x 52
ESPRIT 2V 60/2 25L	25	X	2	210 - 7,4	75 - 2,6 (5 BAR) 50 - 1,8 (8 BAR)	8 - 114	0,75 - 0,5	Single	2800	62	23 - 51	41 x 45 x 52

With automatic thermal protection

Con protezione termica automatica

EXTREME 2D 25 Liters



Packing Imballo

Weight with box (EXTREME 2D 1,5HP) :
44 Kg - 97 Lbs
Box dimensions cm : 58,5 x 57,5 x 66 (0,20m³)
Quantity per (110x110) PALLET : 8 pc.



EXTREME 2V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME 2V 25L 1,50	25	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	33 - 73	46 x 43 x 66
EXTREME 2D 25L 1,50	25	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	41 - 90,5	58,5 x 43 x 66
EXTREME 2V 25L 2,00	25	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	35,5 - 78	46 x 43 x 66
EXTREME 2D 25L 2,00	25	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	44 - 97	58,5 x 43 x 66

With automatic thermal protection
Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 05

Packing

Imballo

Weight with box (EXTREME 2D 1,5HP) :

51 Kg - 112,4 Lbs

Box dimensions cm : 60 x 45 x 86 (0,23 m³)

Quantity per (110x110) PALLET : 8 pc.

EXTREME 2V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50 CV	1 - 2	1	1
EXTREME 2,00 CV	2 - 3	2	1
EXTREME 2,50 CV	3 - 4	3	2

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME 2V 50L 1,50	50	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	39,5 - 87	46 x 43 x 83
EXTREME 2D 50L 1,50	50	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	48 - 106	58,5 x 43 x 83
EXTREME 2V 50L 2,00	50	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	42,5 - 94	46 x 43 x 83
EXTREME 2D 50L 2,00	50	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	51 - 112	58,5 x 43 x 83
EXTREME 2V 50L 2,50	50	X	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,84	Single	73	52 - 115	64 x 45 x 89
EXTREME 2D 50L 2,50	50	✓	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,84	Single	73	60 - 132	64 x 45 x 89

With automatic thermal protection

Con protezione termica automatica

EXTREME 2D 100 Liters



Packing Imballo

Weight with box (EXTREME 2D 2,0HP) :
Kg 70 - 154 Lbs
Box dimensions cm : 60 x 45 x 142 (0,39 m³)
Quantity per (110x110) PALLET : 2 pc.

EXTREME 2V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 2,00	2 - 3	2	1
EXTREME 2,50	3 - 4	3	2

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Type	dB	Kg - Lbs	cm (LxWxH)
EXTREME 2V 100L 2,00	100	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	57 - 127	46 x 43 x 140
EXTREME 2D 100L 2,00	100	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	66 - 145	58,5 x 43 x 140
EXTREME 2V 100L 2,50	100	X	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,84	Single	73	67 - 147	64 x 45 x 140
EXTREME 2D 100L 2,50	100	✓	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,84	Single	73	75 - 165	64 x 45 x 140

With automatic thermal protection

Con protezione termica automatica


MADE IN ITALY BY NARDI COMPRESSORI

Page 07

Packing Imballo

Weight with box (EXTREME 3D 1,5HP) :
50 Kg - 112 Lbs
Box dimensions cm : 88 x 48 x 72 (0,30 m³)
Quantity per (110x110) PALLET : 8 pc.

EXTREME 3V

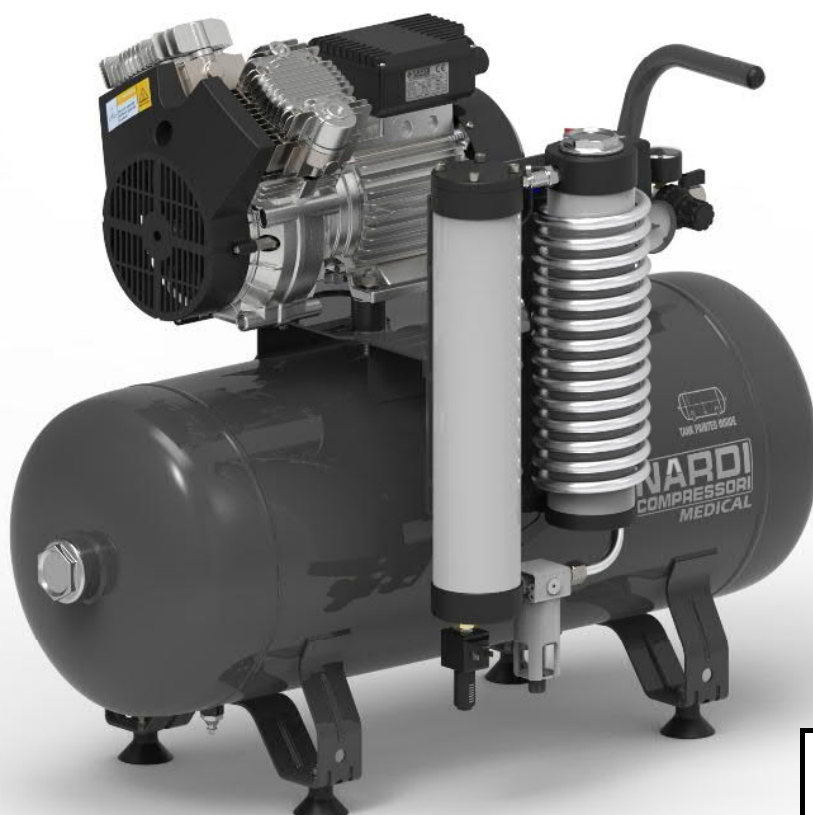

 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME 3V 50L 1,50	50	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	40 - 88	84 x 31 x 68
EXTREME 3D 50L 1,50	50	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	72	47,5 - 105	84 x 44 x 68
EXTREME 3V 50L 2,00	50	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	43 - 95	84 x 31 x 68
EXTREME 3D 50L 2,00	50	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	50,5 - 112	84 x 44 x 68

With automatic thermal protection

Con protezione termica automatica

EXTREME 4D 90 Liters



Packing
Imballo

Weight with box (EXTREME 4D 2,5HP) :
96 Kg - 211 Lbs
Box dimensions cm : 113 x 60 x 84 (0,54m³)
Quantity per EUR PALLET : 2 pc.

EXTREME 4V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 2,00	2 - 3	2	1
EXTREME 2,50	3 - 4	3	2

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME 4V 90L 2,00	90	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	67 - 147	110 x 54 x 81
EXTREME 4D 90L 2,00	90	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	73	76 - 167	110 x 59 x 81
EXTREME 4V 90L 2,50	90	X	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,85	Single	73	77 - 170	110 x 54 x 81
EXTREME 4D 90L 2,50	90	✓	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,85	Single	73	86 - 189	110 x 59 x 81

With automatic thermal protection
Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 09

Packing Imballo

Weight with box (EXTREME 2D SIL 25L 1,5HP) :
45,5 Kg - 100 Lbs
Box dimensions cm : 70 x 43 x 78 (0,23m³)
Quantity per EUR PALLET : 4 pc.



EXTREME 2V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Type	dB	Kg - Lbs	cm (LxWxH)
EXTREME 2V SIL 25L 1,50	25	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	66	38,5 - 85	58 x 43 x 78
EXTREME 2D SIL 25L 1,50	25	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	66	45,5-100	70 x 43 x 78
EXTREME 2V SIL 25L 2,00	25	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	67	41 - 90	58 x 43 x 78
EXTREME 2D SIL 25L 2,00	25	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	67	48 - 106	70 x 43 x 78

With automatic thermal protection

Con protezione termica automatica

EXTREME 2D SIL 50 Liters



Packing Imballo

Weight with box (EXTREME 2D SIL 50L 1,5HP) :
48 Kg - 106 Lbs
Box dimensions cm : 70 x 43 x 98 (0,30m³)
Quantity per EUR PALLET : 4 pc.

EXTREME 2V



Dental unit 	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Type	dB	Kg - Lbs	cm (LxWxH)
EXTREME 2V SIL 50L 1,50	50	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	66	45 - 99	58 x 43 x 98
EXTREME 2D SIL 50L 1,50	50	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	66	53 - 117	70 x 43 x 98
EXTREME 2V SIL 50L 2,00	50	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	67	47 - 103	58 x 43 x 98
EXTREME 2D SIL 50L 2,00	50	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	67	55 - 122	70 x 43 x 98

With automatic thermal protection

Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 11

Packing Imballo

Weight with box (EXTREME SD 1,5HP) :
78 Kg - 172 Lbs
Box dimensions cm : 73 x 47 x 79 (0,24m³)
Quantity per EUR PALLET : 4 pc.



Dental unit 	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Type	dB	Kg - Lbs	cm (LxWxH)
EXTREME SV 30L 1,50	30	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	61	68 - 150	65 x 44 x 76
EXTREME SD 30L 1,50	30	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	61	75 - 165	71 x 44 x 76
EXTREME SV 30L 2,00	30	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	62	71,5 - 158	65 x 44 x 76
EXTREME SD 30L 2,00	30	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	62	78,5 - 173	71 x 44 x 76

With automatic thermal protection
Con protezione termica automatica

EXTREME SD 50 Liters



Packing

Imballo

Weight with box (EXTREME SD 1,5HP) :

90 Kg - 199 Lbs

Box dimensions cm : 85 x 52 x 84 (0,34m³)

Quantity per EUR PALLET : 2 pc.

 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 1,50	1 - 2	1	1
EXTREME 2,00	2 - 3	2	1

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Type	dB	Kg - Lbs	cm (LxWxH)
EXTREME SV 50L 1,50	50	X	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	59	83 - 183	72 x 50 x 82
EXTREME SD 50L 1,50	50	✓	2	240 - 8,5	125 - 4,4 (5 BAR) 100 - 3,5 (8 BAR)	8 - 114	1,50 - 1,1	Single	59	90 - 199	85 x 50 x 82
EXTREME SV 50L 2,00	50	X	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	60	86,5 - 195	72 x 50 x 82
EXTREME SD 50L 2,00	50	✓	2	350 - 12,4	156 - 5,5 (5 BAR) 127 - 4,5 (8 BAR)	8 - 114	2,00 - 1,5	Single	60	95,5 - 210	85 x 50 x 82

With automatic thermal protection

Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 13



Packing
Imballo

Weight with box (EXTREME SD 2,5HP) :
161 Kg - 355 Lbs
Box dimensions cm : 96 x 56 x 110 (0,60 m³)
Quantity per EUR PALLET : 1 pc.

 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 2,50	3 - 4	3	2

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME SV 100L 2,50	100	X	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,85	Single	62	154 - 340	96 x 56 x 102
EXTREME SD 100L 2,50	100	✓	4	480 - 17,0	240 - 8,5 (5 BAR) 190 - 6,7 (8 BAR)	8 - 114	2,50 - 1,85	Single	62	161 - 355	96 x 56 x 110

With automatic thermal protection
Con protezione termica automatica

Line for CAD-CAM use



SPECIFIC LINE FOR CAD-CAM.

(ENG) This line is specifically designed to meet the requirements of this sector, namely huge air demands with a high degree of purity. For this specific application, NARDI COMPRESSORI has developed a special absorption dryer for continuous use, upgraded the size of the motor and improved its ventilation, in order to offer machines that can ensure constant 8 bar and even very intense work cycles.



LINEA DEDICATA PER CAD-CAM.

(IT) Questa linea è appositamente progettata per soddisfare le esigenze di questo settore, ovvero richieste d'aria elevate con un alto grado di purezza. Per questa specifica applicazione NARDI COMPRESSORI ha sviluppato uno speciale essiccatore ad assorbimento per uso continuo, ha rivisto le dimensioni del motore e ne ha migliorato la ventilazione, in modo da offrire macchine che possano assicurare gli 8 bar costanti e cicli di lavoro anche molto intensi.

With automatic thermal protection

Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 15



Packing
Imballo

Weight with box (EXTREME SD 3,0HP) :
77 Kg - 170 Lbs
Box dimensions cm : 74 x 45 x 140 (0,47m³)
Quantity per EUR PALLET : 4 pc.

EXTREME 4V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 3,00	4 - 5	4	3

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tipe	dB	Kg - Lbs	cm (LxWxH)
EXTREME SV 100L 3,00	100	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	73	69 - 152	64 x 45 x 140
EXTREME SD 100L 3,00	100	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	73	77 - 170	74 x 45 x 140
EXTREME SV 100L 3,00	100	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	73	69 - 152	64 x 45 x 140
EXTREME SD 100L 3,00	100	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	73	77 - 170	74 x 45 x 140

With automatic thermal protection

Con protezione termica automatica

EXTREME 4D 90 Liters



Packing
Imballo

Weight with box (EXTREME 4D 3,0HP) :
98 Kg - 211 Lbs
Box dimensions cm : 113 x 70 x 84 (0,58m³)
Quantity per EUR PALLET : 2 pc.

EXTREME 4V



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 3,00	4 - 5	4	3

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME 4V 90L 3,00	90	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	73	67 - 147	110 x 54 x 81
EXTREME 4D 90L 3,00	90	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	73	76 - 167	110 x 59 x 81
EXTREME 4V 90L 3,00	90	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	73	77 - 170	110 x 54 x 81
EXTREME 4D 90L 3,00	90	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	73	86 - 189	110 x 59 x 81

With automatic thermal protection
Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI

Page 17



Packing
Imballo

Weight with box (EXTREME SD 3,0HP) :
167 Kg - 368 Lbs
Box dimensions cm : 96 x 56 x 110 (0,60 m³)
Quantity per EUR PALLET : 1 pc.

Dental unit 	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 3,00	4 - 5	3	2

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME SV 100L 3,00	100	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	62	155 - 342	96 x 56 x 102
EXTREME SD 100L 3,00	100	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Single	62	167 - 369	96 x 56 x 110
EXTREME SV 100L 3,00	100	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	62	155 - 342	96 x 56 x 102
EXTREME SD 100L 3,00	100	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8,8 (8 BAR)	10 - 145	3,00 - 2,20	Three	62	167 - 369	96 x 56 x 110

With automatic thermal protection


Con protezione termica automatica

MULTI-POWER LINE



THE NEW NARDI LP CONTROL DEVICE - TOTAL CONTROL, COMPLETELY INDEPENDENT, ENERGY SAVING EFFICIENCY AND SAFETY RELIABILITY OF NARDI COMPRESSORI.

(ENG) This is an automatic electronic controller, flexible and reliable with a large display and touch commands in an easy-to-use modern and compact design. The LP CONTROL DEVICE controls up to three pump units autonomously depending on the air quantity requirement of the user which rules the starting up of a single or more pump units. Therefore, the amount of air needed would give the system the input to run the pump units one by one to cater for the air quantity needed. When a limited amount of air is needed, the LP CONTROL DEVICE starts a single pump unit, on rotation, this way all pump units will have an equal working time. If more air is needed and a single pump unit is not enough, then the LPCD will start the second electric motor and in case more power is needed to reach the pressure needed, it will also run the third electric motor. This kind of working sequence of the pump units is an important energy saver since it avoids frequent and high peaks of electrical consumption. In

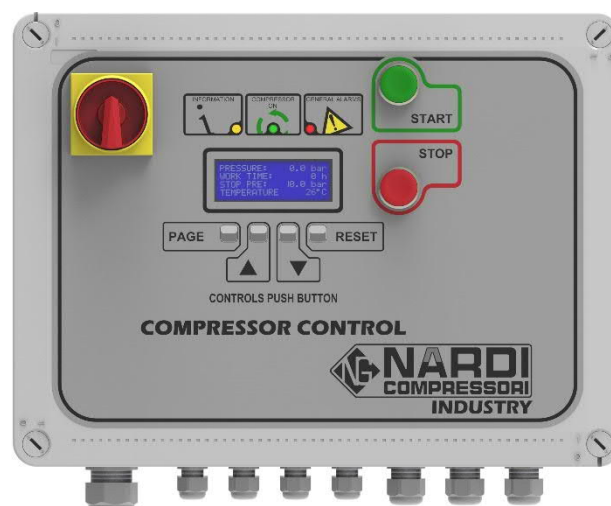

manual mode, during maintenance operations on one or more pump units, the user selects which pump unit to run and which one/s to keep turned off. The automatic control device memorizes the total working hours of the compressor and the working hours of every single pump unit and dryer installed. The motor switches have a thermal control which, in case of an electric overload, stops the motor and informs the automatic system of an ongoing malfunction. This will immediately start another motor this way there is no halt in the work of the user. The LPCD launches a message when the compressor is due for regular maintenance since the working hours of all pump units on the compressor are equal, and also in case any alarm activates. It is possible to collocate up to 10 LP CONTROL DEVICES, therefore a maximum of 30 pump units by means of an RS485 connection. This controlled system consents a better stability of the process, optimization of the compressor and the pumping units, maximization of the working hours, reduction of energy costs and minimum maintenance operations. Upon request, it is possible to add an extra module to have the LPCD send notifications to a determined email address giving real time information about the conditions of the compressor (IoT system).

IL NUOVO NARDI LP CONTROL DEVICE - TOTALE CONTROLLO, COMPLETA AUTONOMIA, EFFICIENTE RISPARMIO ENERGETICO E L'AFFIDABILE SICUREZZA NARDI COMPRESSORI.

(ITA) Sistema di gestione elettronico automatico, flessibile ed affidabile, grande display e comoda pulsantiera, di facile utilizzo e dal design moderno e compatto. Il sistema di gestione LP CONTROL DEVICE comanda fino a tre gruppi pompanti in maniera autonoma, in base alla richiesta d'aria da parte dell'utente. L'avviamento di uno o più gruppi pompanti, è dato dal consumo d'aria dell'utente. In base alla quantità di aria richiesta, la gestione automatica LP CONTROL DEVICE accenderà uno alla volta i motori elettrici in base al fabbisogno d'aria. Se il consumo d'aria è ridotto, farà partire un solo motore, a rotazione con gli altri due, in base al monte ore di lavoro accumulato. Se invece il fabbisogno d'aria è più elevato ed un solo pompante non è sufficiente, allora la gestione automatica farà partire anche il secondo motore elettrico, e di conseguenza, a valori di pressione stabili, accenderà anche il terzo motore elettrico. Questo tipo di gestione dei pompanti, fa risparmiare corrente elettrica ed evita spunti di consumo elettrico frequenti e troppo elevati. In modalità manuale invece, è l'utente a decidere quale gruppo pompante utilizzare e quale invece fermare, in caso di manutenzione ad uno o più dei pompanti.

Il sistema di gestione automatico LP CONTROL DEVICE controlla le ore di lavoro totali del compressore e le ore di lavoro di ogni singolo gruppo pompante e, di conseguenza, anche degli essiccatori. Il comando motori è provvisto di controllo termico, che in caso di eccessivo assorbimento di corrente elettrica, ferma il motore, comunica alla gestione automatica il malfunzionamento, e fa partire subito un altro motore elettrico, così da non fermare l'utente nel suo lavoro. Il sistema avvisa quando il compressore deve fare la manutenzione programmata ed ogni eventuale allarme riscontrato. La manutenzione viene effettuata su tutti i gruppi pompanti durante la stessa manutenzione, in quanto le ore di lavoro sono uguali per tutti i pompanti presenti sul compressore. È possibile affiancare e collegare tra loro fino a 10 LP CONTROL DEVICE, quindi fino a 30 gruppi pompanti, tramite collegamento RS485.

Questo sistema di gestione controllata, consente maggiore stabilità del processo, l'ottimizzazione del compressore e dei gruppi pompanti, l'ottimizzazione dei tempi di lavoro, la riduzione del costo energetico e manutenzione minima. A richiesta, è possibile aggiungere il modulo per la comunicazione dei guasti tramite mail, ad un determinato indirizzo di posta elettronica, per essere sempre informati in tempo reale sulle condizioni del compressore (sistema IoT).

Packing

Imballo

Weight with box (EXTREME MPD 270 9,0HP) :

280 Kg - 617 Lbs

Box dimensions cm : 170 x 95 x 110 (1,78 m³)

Quantity per PALLET : 1 pc.

EXTREME MPV



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 3,00	6	4	3
EXTREME 6,00	12	8	6
EXTREME 9,00	18	12	9

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tip	dB	Kg - Lbs	cm (LxWxH)
EXTREME MPV 270L 3,0	270	X	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8.5 (8 BAR)	10 - 145	3,00 - 2,2	Three	73	116 - 255	160 x 64 x 87
EXTREME MPD 270L 3,0	270	✓	4	600 - 21,2	312 - 11 (5 BAR) 250 - 8.5 (8 BAR)	10 - 145	3,00 - 2,2	Three + N	73	130 - 287	160 x 93 x 87
EXTREME MPV 270L 6,0	270	X	8	1200 - 42,4	324 - 22 (5 BAR) 500 - 18 (8 BAR)	10 - 145	6,00 - 4,4	Three	74	165 - 364	160 x 64 x 87
EXTREME MPD 270L 6,0	270	✓	8	1200 - 42,4	324 - 22 (5 BAR) 500 - 18 (8 BAR)	10 - 145	6,00 - 4,4	Three	74	193 - 425	160 x 93 x 87
EXTREME MPV 270L 9,0	270	X	12	1800 - 63,6	936 - 33 (5 BAR) 750 - 27 (8 BAR)	10 - 145	9,00 - 6,6	Three	75	198 - 436	160 x 64 x 87
EXTREME MPD 270L 9,0	270	✓	12	1800 - 63,6	936 - 33 (5 BAR) 750 - 27 (8 BAR)	10 - 145	9,00 - 6,6	Three	75	240 - 530	160 x 93 x 87

With automatic thermal protection

Con protezione termica automatica

EXTREME MP 500 Liters



Packing

Imballo

Weight with box (EXTREME MPD 500 9,0HP) :

340 Kg - 750 Lbs

Box dimensions cm : 200 x 120 x 120 (2,9 m³)

Quantity per PALLET : 1 pc.

EXTREME MPV



 Dental unit	Pneumatic		
	Electric-Pneumatic		
	Electric		
EXTREME 6,00	12	8	6
EXTREME 9,00	18	12	9

Compressor Model	Tank	Dryer	Cylinders	Air Intake	Air delivery	Pressure	Power	Phase	Noise	Weight	Dimensions
Modello Compressore	Serb.	Essicat.	Cilindri	Aria Aspirata	Aria resa	Pressione	Motore	Fase	Rumore	Peso	Dimensioni
Type	Liters	Complete	N°	l/min. - C.F.M.	l/min. - C.F.M.	Bar - Psi	HP - kW	Tipe	dB	Kg - Lbs	cm (LxWxH)
EXTREME MPV 500L 6,00	500	X	8	1200 - 42,4	324 - 22 (5 BAR) 500 - 18 (8 BAR)	10 - 145	6,00 - 4,4	Three	74	220 - 485	195 x 64 x 110
EXTREME MPD 500L 6,00	500	✓	8	1200 - 42,4	324 - 22 (5 BAR) 500 - 18 (8 BAR)	10 - 145	6,00 - 4,4	Three	74	248 - 547	195 x 93 x 110
EXTREME MPV 500L 9,00	500	X	12	1800 - 63,6	936 - 33 (5 BAR) 750 - 27 (8 BAR)	10 - 145	9,00 - 6,6	Three	75	253 - 558	195 x 64 x 110
EXTREME MPD 500L 9,00	500	✓	12	1800 - 63,6	936 - 33 (5 BAR) 750 - 27 (8 BAR)	10 - 145	9,00 - 6,6	Three	75	295 - 650	195 x 93 x 110

With automatic thermal protection

Con protezione termica automatica

MADE IN ITALY BY NARDI COMPRESSORI


Page 21

DRYER





Model	Air filtration (5 BAR)	Duty cycle	Power	Purge	Dew Point	Max Pressure	Weight	Dimensions
Modello	Aria filtrata (5 BAR)	Ciclo di lavoro	Consumo	Perdita di carico	Dew Point	Pressione mass.	Peso	Dimensioni
Type	l/min. - C.F.M. - m³/h	%	Watt	%	°C	Bar - Psi	Kg - Lbs	cm (LxWxH)
DRYER 400	400 - 14 - 24	80 % (20 min)	18	1 %	-20	10 - 145	7 - 15,5	51 x 26 x 13
DRYER D600	600 - 21,2 - 36	100 %	46	6 %	-40	10 - 145	12 - 26,5	50 x 33 x 26
DRYER M400	400 - 14 - 24	100 %	10	14 %	-5	10 - 145	7 - 15,5	57 x 27 x 26

FILTER DEPURATOR

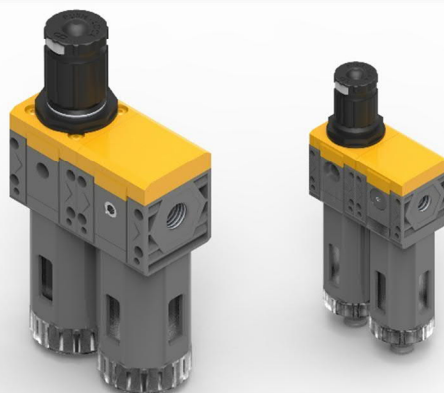


Model	Air filtration (8 BAR)	Duty cycle	Connection	Depuration	Max Pressure	Weight	Dimensions
Modello	Aria filtrata (8BAR)	Ciclo di lavoro	Attacco	Purificazione	Pressione mass.	Peso	Dimensioni
Type	l/min. - C.F.M. - m³/h	%	in / out	µm	Bar - Psi	Kg - Lbs	cm (LxWxH)
DEPURATOR 1/4"	460 - 14 - 24	100 %	G1/4" - G1/4"	5 + 0,01	10 - 145	2 - 4,5	8,5 x 4,5 x 19
DEPURATOR 1/2"	990 - 14 - 24	100 %	G1/2" - G1/2"	5 + 0,01	10 - 145	5 - 11	12 x 6 x 23,5

PUMP UNIT

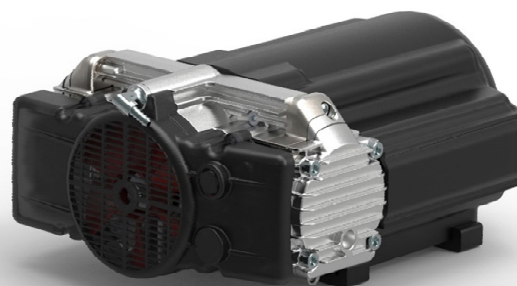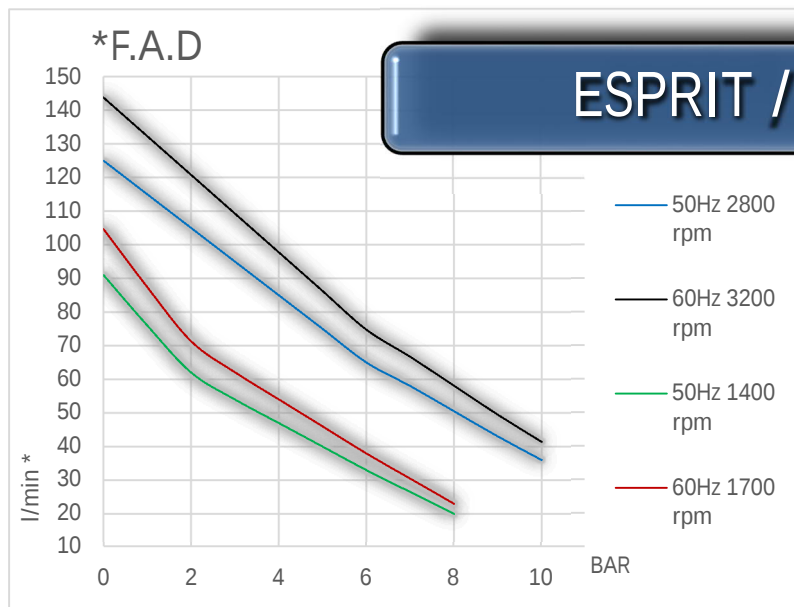

Compressor Model	Cylinders	Air Intake	Pressure	Power	Phase	Revolution	Noise	Weight	Dimensions
Modello Compressore	Cilindri	Aria Aspirata	Pressione	Motore	Fase	Giri	Rumore	Peso	Dimensioni
Type	N°	l/min. - C.F.M. - m³/h	Bar - Psi	HP - kW	Tip	Rpm	dB	Kg - Lbs	cm (LxWxH)
ESPRIT UNIT 60/4	2	105 - 3,7 - 6,3 (50 Hz) 125 - 4,4 - 7,5 (60 Hz)	8 - 114	0,75 - 0,5	Single	1400 1700	48	8 - 17	29 x 22 x 14
ESPRIT UNIT 60/2	2	210 - 7,4 - 12,6 (50 Hz) 250 - 8,8 - 15 (60 Hz)	8 - 114	0,75 - 0,5	Single	2800 3300	55	8 - 17	29 x 22 x 14

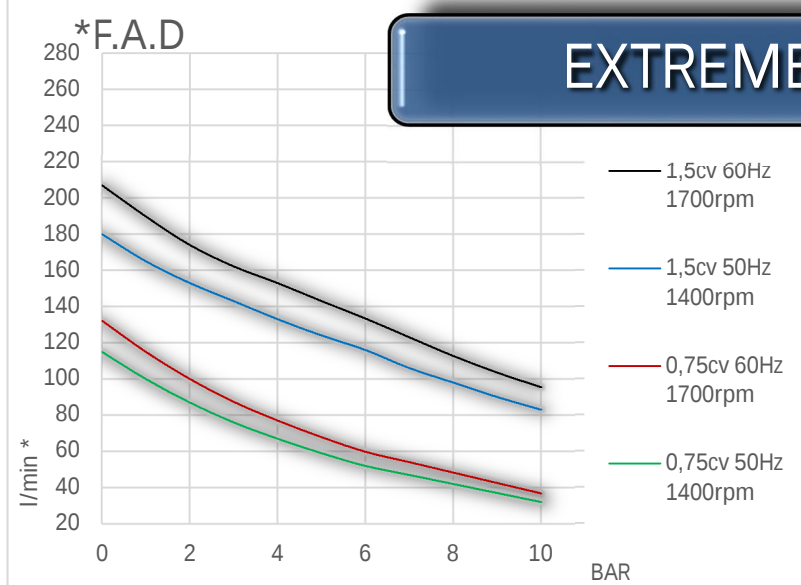


Compressor Model	Cylinders	Air Intake	Max Pressure	Power	Revolution	Noise	Weight	Dimensions
Modello Compressore	Cilindri	Aria Aspirata	Pressione mass.	Motore	Giri	Rumore	Peso	Dimensioni
Type	N°	l/min. - C.F.M. - m³/h	Bar - Psi	HP - kW	Rpm	dB	Kg - Lbs	cm (LxWxH)
EXTREME UNIT 0,75 CV	2	200 - 7,0 - 12 (50 Hz) 240 - 8,5 - 14,4 (60 Hz)	10 - 145	0,75 - 0,55	1400 1700	72	14 - 31	39 x 28 x 26
EXTREME UNIT 1,50 CV	2	240 - 8,5 - 14,4 (50 Hz) 290 - 10,1 - 17,2 (60 Hz)	10 - 145	1,50 - 1,1	1400 1700	72	16,5 - 36	39 x 28 x 26

With automatic thermal protection - * F.A.D. values have a tolerance $\pm 5\%$


Con protezione termica automatica - * I valori indicati hanno una tolleranza $\pm 5\%$

MADE IN ITALY BY NARDI COMPRESSORI

Page 23

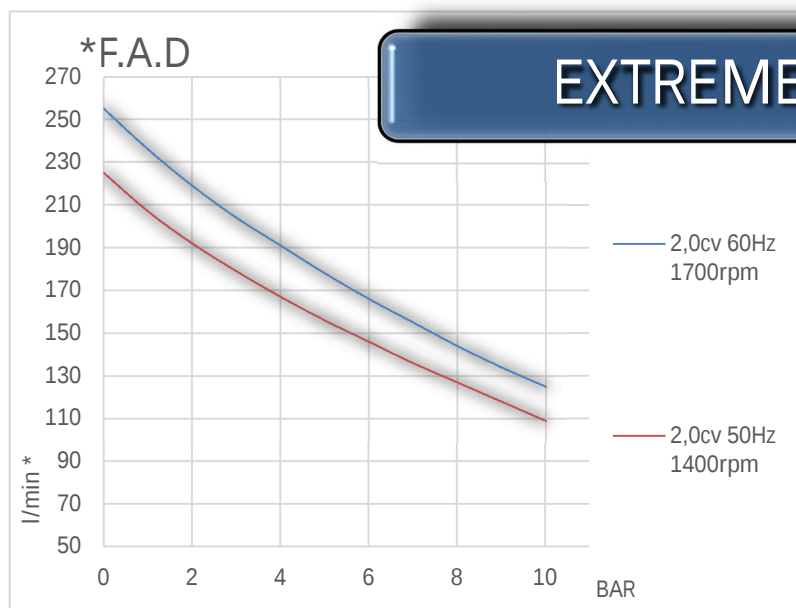

Compressor Model	Cylinders	Air Intake	Max Pressure	Power	Revolution	Noise	Weight	Dimensions
Modello Compressore	Cilindri	Aria Aspirata	Pressione mass.	Motore	Giri	Rumore	Peso	Dimensioni
Type	N°	l/min. - C.F.M. - m³/h	Bar - Psi	HP - kW	Rpm	dB	Kg - Lbs	cm (LxWxH)
EXTREME UNIT 2 CV 1400	2	350 – 12,4 – 21 (50 Hz) 420 – 14,8 – 25 (60 Hz)	10 – 145	2,00 – 1,5	1400 1700	73	20 – 44	43 x 28 x 29

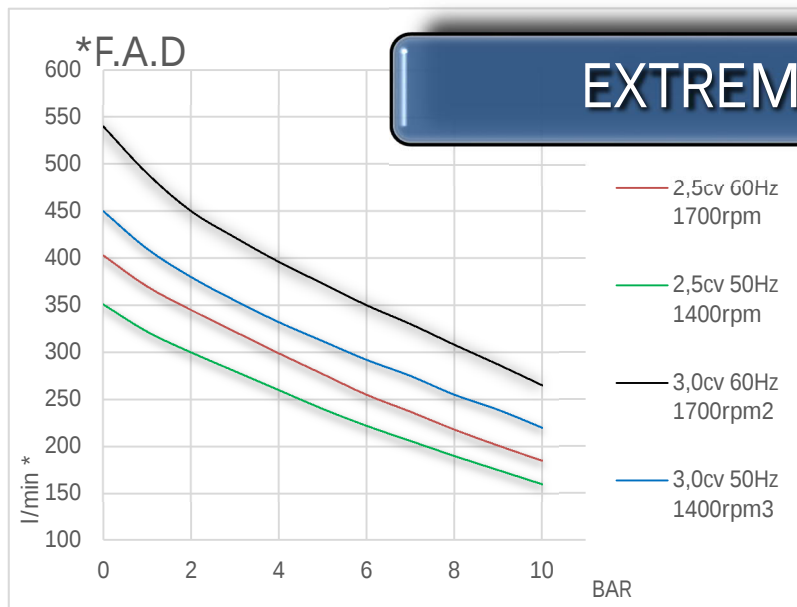

Compressor Model	Cylinders	Air Intake	Max Pressure	Power	Revolution	Noise	Weight	Dimensions
Modello Compressore	Cilindri	Aria Aspirata	Pressione mass.	Motore	Giri	Rumore	Peso	Dimensioni
Type	N°	l/min. - C.F.M. - m³/h	Bar - Psi	HP - kW	Rpm	dB	Kg - Lbs	cm (LxWxH)
EXTREME 4 UNIT 2,50 CV	4	480 – 17 – 29 (50 Hz) 580 – 20 – 35 (60 Hz)	10 – 145	2,50 – 1,84	1400 1700	73	30 – 66	61 x 28 x 31
EXTREME 4 UNIT 3,00 CV	4	600 – 21,2 – 36 (50 Hz) 720 – 25,5 – 43 (60 Hz)	10 – 145	3,00 – 2,2	1400 1700	73	33 – 73	64 x 28 x 31

With automatic thermal protection - * F.A.D. values have a tolerance $\pm 5\%$

Con protezione termica automatica - * I valori indicati hanno una tolleranza $\pm 5\%$

General terms



Standard Delivery

Fornitura standard

Standard Volt/Hertz/Phase : 230/50/1
Other : 100⇒240 Volt/50-60 Hz/1 Phase

STANDARD COLOR : GREY
Other colors min. delivery : 50 pcs

Standard plug :
SCHUKO

Standard fitting :
1/4"G BALL VALVE

General terms of sale (ENGLISH)

- This catalogue cancels and substitutes the previous versions.
- Once order is confirmed all terms of sale are implicitly accepted.
- The goods are always delivered ex works and at the buyer's risk.
- Any possible complaints about the packing or the delivery are accepted only if notified within 8 days from the date of receipt.
- Our production is warranted for 24 months from the delivery of the goods. Warranty is excluded for improper use and consequent damages. Warranty is guaranteed only for manufacturing fault.
- Transport costs are always and completely at customer charge.
- Delivery terms are always to be considered as approximate. Shortage of all or part of the goods on delivery can be claimed only if accompanied by written terms agreed upon. Otherwise redraft cannot be claimed.
- All data and technical aspects cited in this catalogue are not binding and may be altered without notice due to technical or manufacturing necessities.
- The invoice referring to the delivered goods is to be paid according to the terms and in the place agreed. Any delay in paying it will enforce late charges.
- In case of disputes, jurisdiction will be the court of Vicenza – Italy.
- This catalogue belongs exclusively to the firm NARDI COMPRESSORI situated in Montecchio Maggiore (VI) Italy; reproduction or copying of the goods without the firm's permission is against the law.

Condizioni generali di vendita (ITALIANO)

- Il presente catalogo annulla e sostituisce il precedente.
- Con il conferimento dell'ordine si considerano accettate tutte le condizioni di vendita.
- La merce si considera franco nostra fabbrica. La merce viaggia sempre a rischio e pericolo dell'acquirente che ha l'obbligo di ritirarla.
- Non si accettano reclami su imballi o spedizioni se non fatti entro 8 giorni dalla consegna della merce.
- I nostri prodotti sono coperti di garanzia per 24 mesi dalla consegna della merce. Non sono garantiti i danni provocati da usura, cattivo uso e non ci assumiamo nessuna responsabilità per danni diretti o indiretti derivanti da tale cattivo uso. La garanzia è valida solo per difetti di fabbricazione.
- Tutti i costi di trasporto sono sempre e completamente a carico del cliente.
- I termini di consegna non sono mai vincolanti. La mancata consegna di parte o tutta l'ordinazione non può dar luogo ad alcuna rivalsa di danni verso la nostra ditta, salvo diversi accordi scritti e regolarmente accettate.
- Tutti i dati e le indicazioni tecniche citate nel presente catalogo non sono impegnativi e potranno essere modificati in caso di necessità tecniche o produttive, senza alcun preavviso.
- Per pagamenti effettuati oltre la scadenza pattuita saranno addebitati interessi moratori al tasso corrente bancario.
- In caso di controversia è riconosciuta la sola competenza del tribunale di Vicenza.
- Questo catalogo è di proprietà esclusiva della ditta Nardi Compressori Srl con sede a Montecchio Maggiore (VI) e senza il consenso della stessa ditta non può essere né copiato, né riprodotto anche solo parzialmente sotto pena delle sanzioni di legge in materia.


This catalogue replaced the previous one

Questo catalogo sostituisce i precedenti

MADE IN ITALY BY NARDI COMPRESSORI

Page 25

Nardi Compressori ITALY

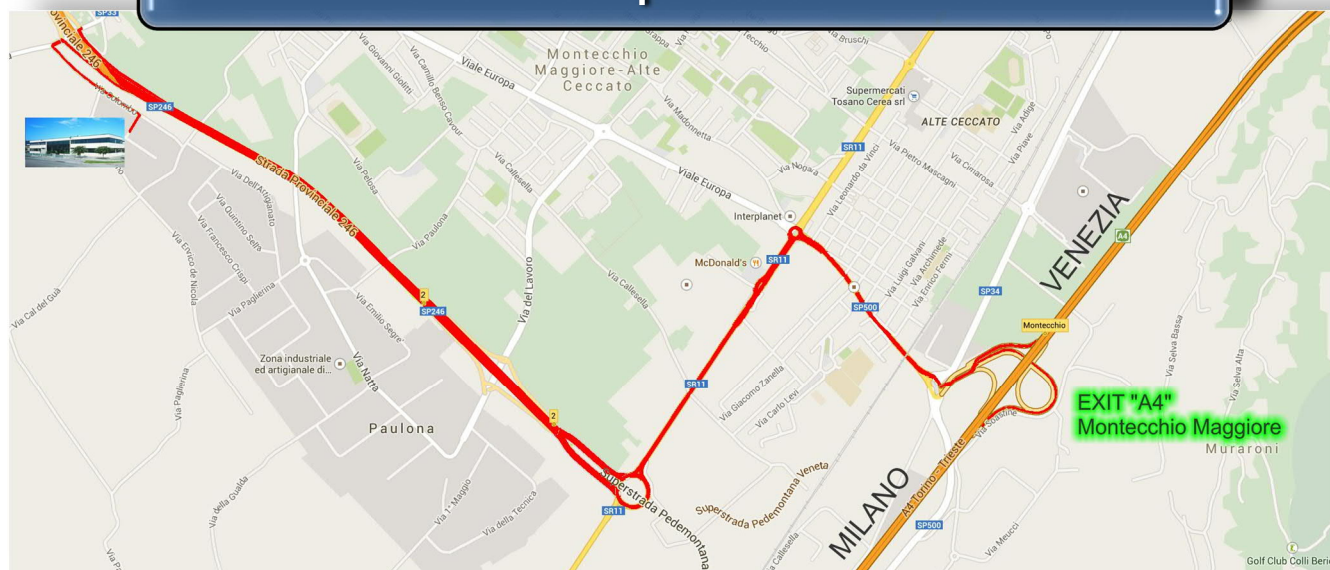

Connection



Railway station VICENZA

15 Km



Airport VERONA "Villafranca"

VRN

55 Km



Airport VENICE "Marco Polo"

VCE

83 Km



Airport MILAN "Orio al Serio"

BGY

150 Km



Airport BOLOGNA "Marconi"

BLQ

170 Km



Airport MILAN "Linate"

LIN

200 Km




Airport MILAN "Malpensa"

MLX

230 Km

PRODUCTION

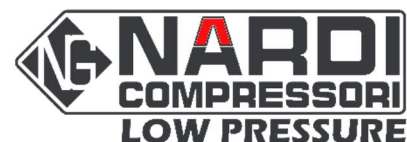

The OIL-LESS series offers a selection of single-phase and three-phase piston compressors. The professional compressors of this line are reliable machines, built with the best materials, to guarantee our customers long life durability with very little maintenance. These compressors are available in 8-10Bar, 15 bar, 21bar, up to a maximum of 30bar and they can be ordered with different frames and tank sizes to satisfy the market demand.



The NARDI Breathing air line has a wide selection of compressors up to 450 Bar for the filling of breathing air cylinders for different uses like: commercial and recreational diving, fire and rescue operations and military applications.



The INDUSTRY line offers a vast selection of compressors suitable for applications in different sectors and industrial markets. All the machines in this range are built with the best materials, to offer robust and durable solutions that can withstand also intense industrial work cycles.



High pressure compressors for refueling of vehicles that run on natural gas. Compressors in this line have been designed for industrial and domestic use, creating simple and economic small refueling stations for private use using gas supplied through the network lines.





Nardi Compressori Srl
Via Marco Polo, 2 36075
Montecchio Maggiore (VICENZA) ITALY
Tel. +39.0444.159111 - 497514 -694885
Fax. +39.0444.159122
Mail. info@nardicompressori.com



www.nardicompressori.com



www.facebook.com/nardicompressorisrl/



www.linkedin.com/company/nardi-compressori-srl/mycompany/

MADE IN VENETO



EXCELLENCE & HIGH TECHNOLOGY